



rauberger
vetmedical



TT-1502 LCD

9.7" LCD Digital Microscope,
5.0M, Android Pad

9.7" LCD / Android Pad / 5.0M

Explore Micro World Together

Multipul Usage / Easy To Share / Nice For Teaching



Excellent Optical System

Fully Coated Lens
Large Eyepiece
Infinity Optional



TRINOCULAR

Trinocular



LED

LED



TRANSMIT

Transmit Light



KOHLER

Kohler Illumination



USB

USB



LCD

LCD



ABBE


Abbe Condenser



3
YEARS

3 YEARS

3 Years Warranty

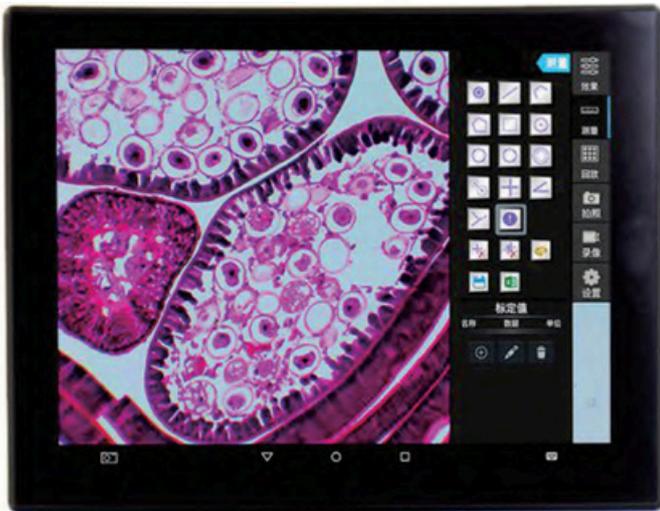



9.7" Android Pad + 5.0M Camera

Sharp Image | All-in-One | Photo & Video




Measuring Software Included

Touch Screen | Full Range Measure | 90 Languages

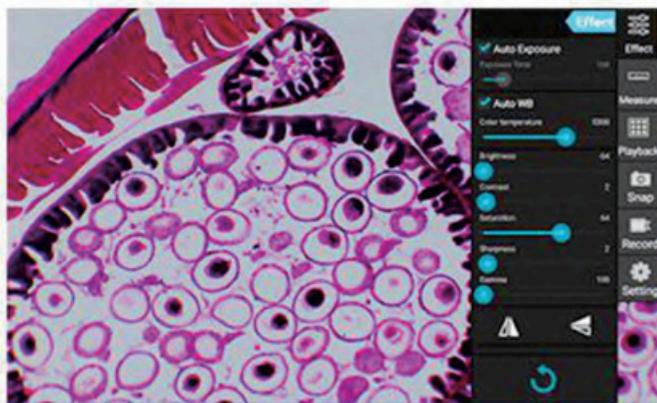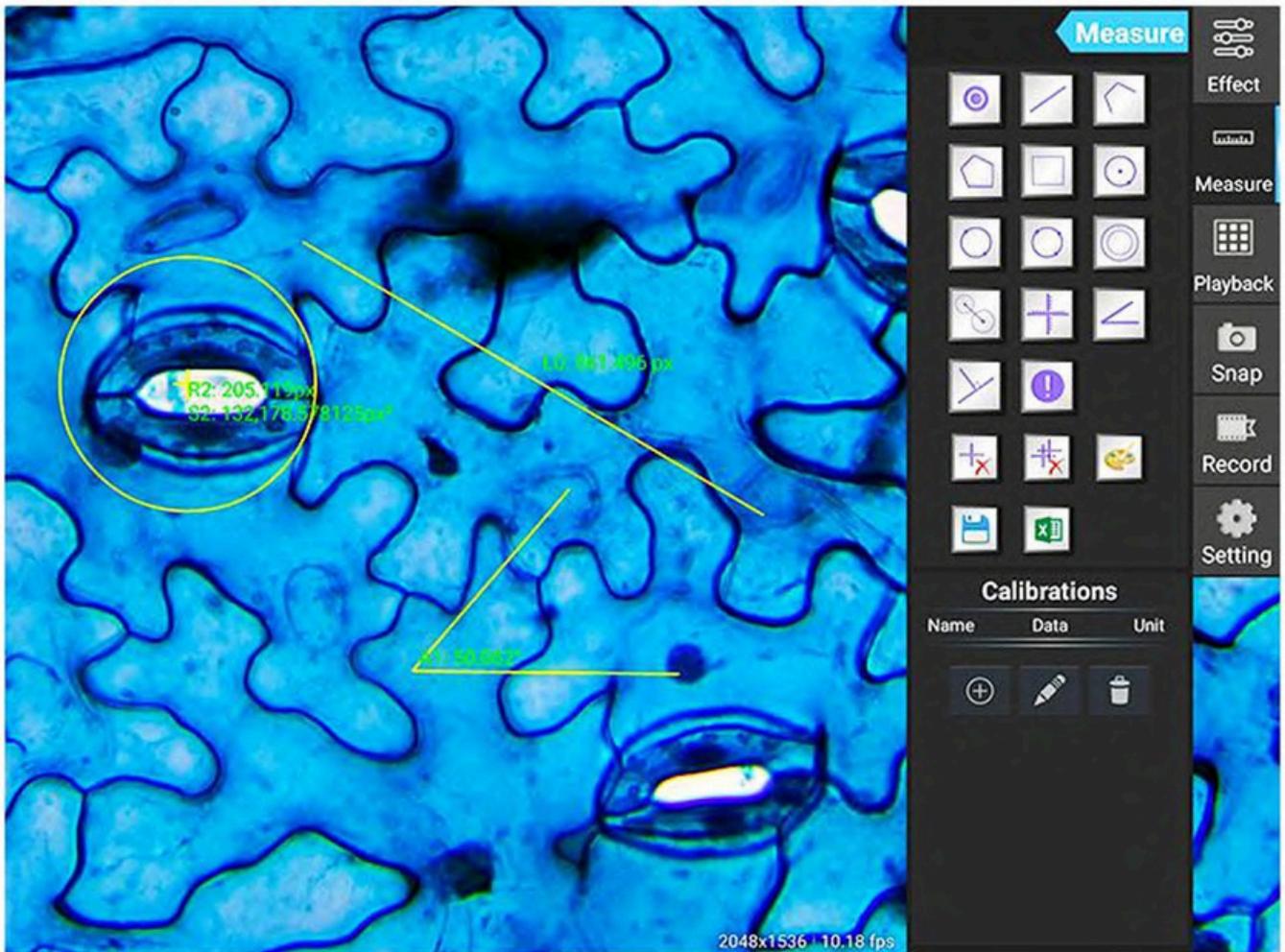

All Optical Parameter Adjustable!

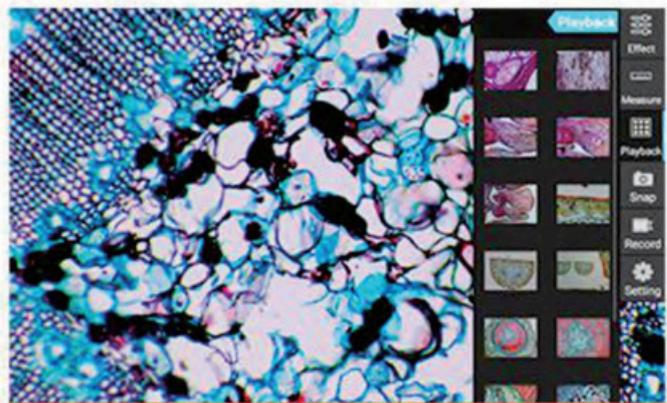

Capture Photo & 1080p Video!



USB / HDMI / TF Card

Product Details »»



High Point View & Wide Field Eyepieces

WF10X/20mm (A12.1502-AB/AT);
WF10X/22mm (A12.1502-B/T)
Diopter Adjustment Range D-8 to +5.
Easy To Find Objectives By Wide Field Vision.



Rotatable Binocular Viewing Head


360 Degree Rotatable Binocular Head.
Interpapillary Distance Adjustment:48mm-75mm.



Coaxial Coarse and Fine Adjustment, With Safety Stop Screw, Tension Adjustable



Achromatic Optical System
Large Size Double Layer Mechanical Stage 142x132mm,
Graphite Surface

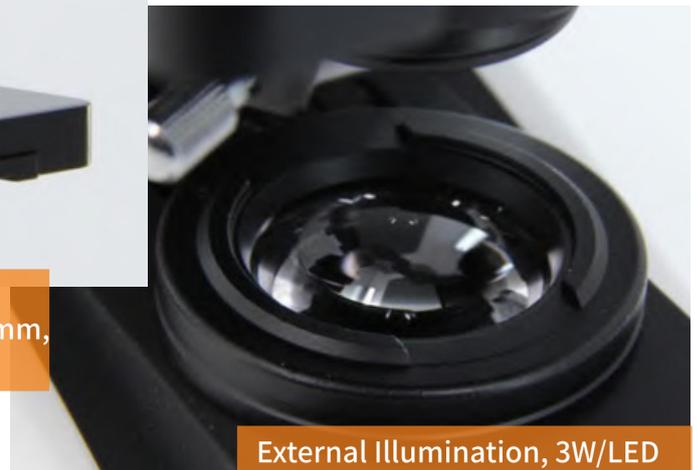

External Illumination, 3W/LED



Abbe Condenser NA1.25, With Iris Diaphragm, Center Adjustable



3W LED Light Source, Removable Lighthouse For Easy Changing Bulb



Switch & Power Socket AC100V-240V
Newly Designed Wire Winding Device For Easily Store Power Cord

Accessories >>>

A51.1502-16



A51.1502-20



A55.1502-05



A5D.1502-D



A5D.1502-I



A5P.1501



Product Specification »»

TT-1502 LCD 9.7" LCD Digital Biological Microscope	
Head	Trinocular Head, 30° Inclined, Interpupillary Distance 48-75mm, Both Eyepiece Tube Diopter Adjustable
Eyepiece	WF10x/20mm, Diopter Adjustable, Dia. 30mm
Objective	High Contrast Chromatic Free 4x, 10x, 40x(S), 100x(Oil, S)
Focus Knob	Coaxial Coarse & Fine Focusing, Fine Division 0.002mm, Coarse Stroke 37.7mm/Rotation, Focusing Range 24mm
Nosepiece	Quadruple, Backward With Click Stops
Stage	Double Layer Mechanical Stage 142x132mm, Graphite Surface, Moving Range 75x50mm
Condenser	Abbe Condenser N.A.1.25, Center Adjustable
Light Source	3W LED Brightness Adjustable
Adapter	C-Mount 1.0x
Package	53*32.5*53.5cm/11.5Kg, 1 Set/Carton
9.7" LCD Pad Digital Camera (A59.3520)	
LCD	9.7" HD IPS Touch Screen With Backlight
Camera	5.0M 1/2.5" CMOS Digital Camera
Resolution	2048x1536 @15FPS
Pixel	2.2umx2.2um
Data Storage	TF Card or Built-in Memory
Data Output	USB2.0, HDMI, TF Card
Photo	5.0M .JPEG
Video	1080p Video Record
System	Android 5.1 System
	90 Languages
	Support Mouse + Keyboard Control
	Support WIFI
Function	Support Bluetooth
	S-Eye Software & Instruction Manual Installed
	Mark: Point, Crosshair, Coordinate, Text
	Length Measure: Straight Line Length, Broken Line Length, Curve Length, Parallel Line Distance, Point To Line Distance
	Distance Between 2 Circle Centers: Radius Fixed Circle, Two-point Fixed Circle, Three-point Fixed Circle
	Geometric Measure: Line Segment Length, Radius Fixed Circle, Two-point Fixed Circle, Three-point Fixed Circle, Concentric Circle
Geometric Area Measure:Polygon, Square	
Power	DC 12V/2A Power Adaptor, AC 100~240V
Size	Pad 238x182x12mm, Camera 40x39mm